

**AFT**  
Flexion Automatique



**Forsco Labs**  
Distributeur matériel - HILL Laboratories

# HILL AFT™ MANUAL DROPS

## Motorized Automatic Flexion

The Hill Automatic Flexion Table is perfect for doctors who want full-featured automatic flexion without the need for manual flexion. The AFT is surprisingly affordable for a table that features variable speed flexion, electric adjustable height and an impressive list of other standard features. Make the AFT even more versatile by adding up to four manual drops and several other options to customize the AFT for your specific needs.



Shown in Navy Upholstery with Beige Base.


Picture with several options

### STANDARDS FEATURES

- Variable speed motorized flexion
- Tilting Headpiece, 30° negative and positive tilt
- Thoracic release and lock
- Pelvic section has lateral flexion and can be locked at any angle up to 21°
- Removable tiller bar at foot-end for lateral movement
- Slide-out ankle support with ankle harnesses
- Adjustable patient gripper bar
- Patient Shut-Off Switch
- Electric height controlled by Dual Rocker Foot Pedals
- Width 24» / Height 21.5» to 29», higher heights available
- Foam Top - 2 1/2» Medium (standard), Soft or Firm
- Choice of 20 Colors for the top cushions
- Beige, Gray or Black Base, enclosed for safety
- One-year warranty on all parts and labor



The AFT has adjustable height 21.5» to 29»

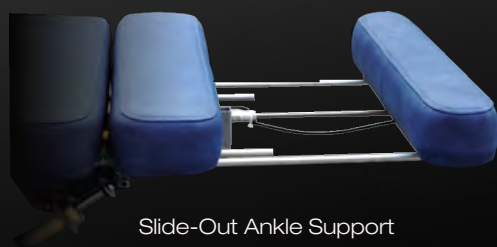

Slide-Out Ankle Support

# AFT MANUAL STANDARDS

All of the following features are standard on the AFT

## ▶ AUTOMATIC MOTORIZED FLEXION

### Flexion variable speed.

The Hill AFT has a Flexion Control Panel that is conveniently located at the foot-end of the table which allows you to set the flexion speed from 1 to 28 rpm (an optional 30 minute timer can also be added). A selector switch allows you to stop the table in the horizontal or flexed position. In the flexed position, with the table at the lowest height, patients can easily enter or exit the table from the foot end. The pelvic section of the AFT moves laterally and can be easily maneuvered during flexion using the Tiller Bar. It can also be locked laterally to either side.

Tilting




Hand lever releases the headpiece to lock at any angle from 30° positive to 30° negative.

Tiller Bar



Inserts for using manual lateral movement or circumduction during a flexion treatment.

Patient Shut-off

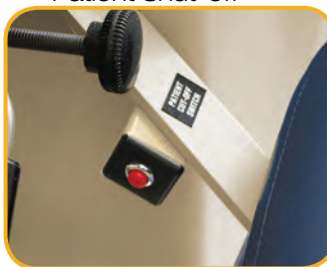

Equipped with a patient shut-off switch located near the left armrest.

Rocker Foot Pedal



Mounted on both sides of the table base for easy height adjustment.



Thoracic release



Patient Gripper Bar



Flexion control

## Main options to complement and customize your AIR-AFT



Flexion and rotation Headpiece



Hill AFT Scoliosis Treatment Package




Minuterie



Raised Headpiece



Gas-powered Stool

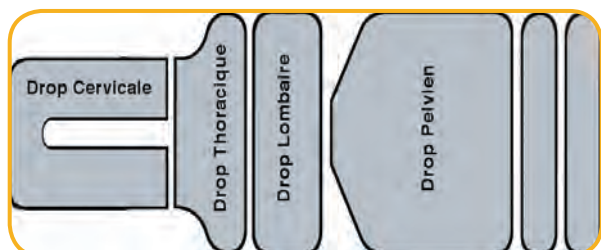

Other complementary options on [www.forscolabs.com](http://www.forscolabs.com)

# SPECIFICATIONS

## AIR-AFT

Electric actuator 60 Hz, 110 volts, 4.8 amps./AirAFT -10 amps. (220-50 cycle also available)

Lifting capacity -450 lbs.

Length - 5' 5» (extends to 6'2» ), with cervical drop - 5'7» to 6'4» • Width - 24»

Range of height; 21 1/2» to 29»; higher heights available

Medium density foam top; firm and soft density upon request

Choice of colors for upholstered top (see below)

Beige, gray or black telescoping base skirting constructed of fire-rated thermo-plastic

Powder-coated internal and external metal parts

Shipping weight 250-300 lbs.


One-year warranty on all parts and labor

## AVAILABLE COLORS

Cushion colors



Base Colors

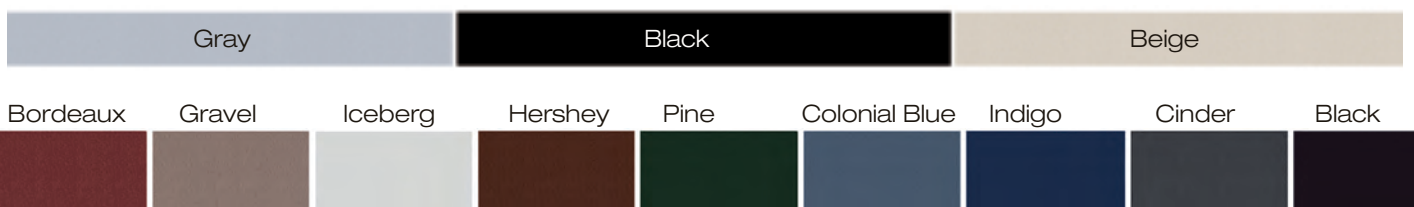

Upholster your table in durable Boltaflex® with Prefixx™, scientifically formulated to give superior protection from staining and abrasions. Finesse and custom colors are available at additional cost. Shades may vary from actual material.

Color choices may change without notice.

Visit [www.forscoLabs.com](http://www.forscoLabs.com) to set your table and make your quote online.

Call one of our staff who will be happy to answer all your questions.

**ForscoLabs**  
Exclusive Distributor - HILL Laboratories

149 avenue du Maine 75014 Paris, France • Tel: +33 145 450 301 • Fax: +33 972 447 924

[contact@forscolabs.com](mailto:contact@forscolabs.com) • [www.forscoLabs.com](http://www.forscoLabs.com)